Orthodontic Elastomeric Separators

Intended Use

Elastomeric separators are placed between adjacent teeth to assist in creating space between them.

Description of the device

Elastomeric separators are typically small, toric-shaped elastomers attached to a handle or in a bulk bag. An alternate product is a barbell shaped solid strip that can be cut to customized length.

Separators are malleable and easily placed between the contact areas of adjacent teeth. They are also referred to as “seps.”

1. Part Numbers

Product	Part	Variants	Material	Composition
Radiopaque Separators	500* 325-*	Blue, Purple	Polyurethane	Thermoplastic Polyurethane

Warnings and Precautions

	Devices are designed and manufactured for single patient use only. The reuse or reprocessing may cause cross infection or loss of mechanical properties due to natural wear of the product.
	Separator breakage, dislodgement or impingement requires immediate attention by the orthodontist.
	Care should be taken to avoid habits or foods that could break or dislodge the separators. Follow orthodontist recommendations for proper at home care.

Clinical Benefits

Elastomeric separators create space between teeth to prepare for placement of orthodontic appliances or to allow visualization of contact areas for recontouring or interproximal reduction (IPR).

Separators, in conjunction with orthodontic appliances, facilitates an aesthetic improvement with the alignment and leveling of the teeth, better hygiene, more efficient masticatory function facilitating swallowing, better digestion, prevention of occlusal problems, better speech and breathing.

Residual Risks

Orthodontic treatments are subject to relapse. Elastomeric separators, as part of the fixed orthodontic system, can contribute to discoloration or decalcification of teeth, bone and root resorption, and periodontal complications.

Potential Side Effects or Adverse Events

- Discomfort and pain associated with separator placement
- Gingival impingement
- Ingestion
- Aspiration
- Elastomeric separators, as part of the fixed orthodontic system, can contribute to difficulties in maintaining oral hygiene, difficulty speaking or chewing and halitosis

Instructions for use

The use of elastomeric separators are restricted to dental professionals with qualification in orthodontics. Use by persons outside this context may result in irreversible damage to the patient.

1. Select the elastomeric separator to be used, or cut customized one from separator strip. Consider radiopacity when selecting a separator.
2. Place the elastomeric separator such that it seats securely at the contact area.
3. Leave in place until ideal separation is achieved. This usually occurs within 3 to 7 days.
4. Patient should be monitored regularly when undergoing orthodontic treatment.

Any serious incident experienced with the use of the product should be reported to the manufacture and the competent authority.

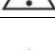
Maintenance and Storage instructions

Store elastomeric separator in product packaging until use.

Disposal

The professional should use good clinical practices for the disposal of medical products in accordance with current sanitary regulations.

Patient Instructions

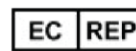
	Caution should be taken not to allow separator to remain in place for longer than 10-14 days.
	Practice good oral hygiene to prevent plaque installation and subsequent periodontal disease.
	Do not floss the areas where separators are in place, dislodgement or impingement could occur.
	Patients should wear a mouthguard when engaging in sports activities; injury to the patient and /or breakage or dislodgment of the appliance can occur. In the event of injury or ingestion, instruct the patient to seek medical care.
	Separator breakage, dislodgement or impingement requires immediate attention by the orthodontist.
	Care should be taken to avoid habits or foods that could break or dislodge the elastomeric separator. Follow orthodontist recommendations for proper at home care.

Warranty

Limited Warranty. Product(s) manufactured by G&H are warranted to be free of defects in materials and manufacture for ninety (90) days from date of delivery. Any material or manufacture defect covered by this limited warranty which occurs during normal use and is reported to Seller in writing during the period of ninety (90) days from the date such products are shipped to customer. G&H's obligation hereunder, upon verification of the defect or error, shall be to provide one of the following, in G&H's sole discretion: (i) replacement at no charge to customer; or (ii) credit the purchase price to customer. If, upon the inspection of G&H product to which this warranty applies, G&H determines that a claimed defect was not due to its manufacture or materials, G&H will invoice customer for the value of product and services at customer's expense. This warranty shall not apply to any product that has been subject to misuse, misapplication, neglect (including but not limited to use of unauthorized parts or attachments), adjustment or repair and/or improper handling and storage. Identifiable items manufactured by others but installed in or affixed to G&H products, or distributed by G&H to customer, are not warranted by G&H but bear only the express warranty, if any, of the manufacturers thereof. No representative of G&H has any authority to waive, alter, vary or add to the terms hereof without prior approval in writing, to customer, signed by an officer of G&H. G&H's liability for its products, whether for breach of contract, negligence, strict liability in tort, or otherwise, shall be limited, in G&H's sole discretion, to the replacement of the products or parts thereof or a credit to the purchase price to customer. G&H will not be liable for any other injury, loss, damage or expense, whether direct or consequential, including but not limited to loss of use, income, profit or production, or increased cost of operation, or spoilage of or damage to material, arising in

connection with the sale, installation, use of, inability to use, or the repair or replacement of, or late delivery of, G&H products.

Limited Warranty. Product(s) manufactured by Highland Metals are warranted to be free from defects. Products must be returned within 90 days from the date of purchase and must be returned in original packaging and in resalable condition. HMI™ reserves the right to final determination as to acceptance of product exchange or credit. Restock fee may apply. Prior to the return of any merchandise, please contact HMI™ Customer Service for a return authorization number

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